




Spectro 20D Plus RS-232

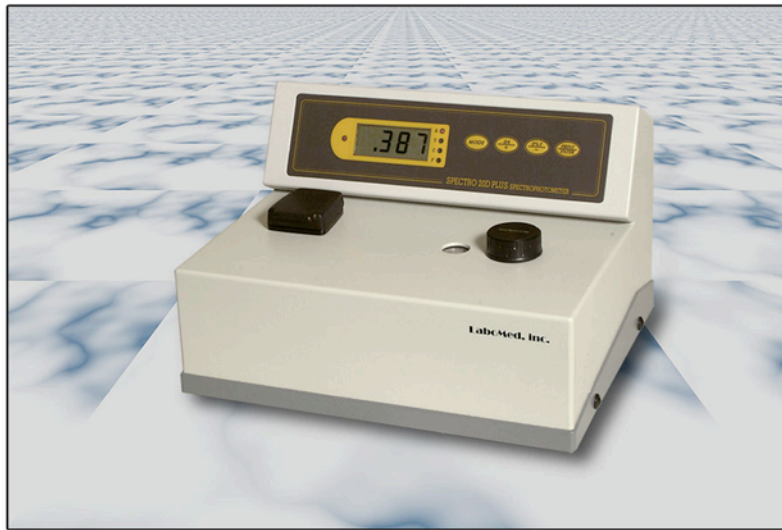

Spectro 20D Plus RS-232C is a traditional analytical device used in small and educational laboratories. This spectrophotometer delivers enhanced ease-of-use, precision and accuracy resulting in time and cost savings, as well as unprecedented confidence in test results. Spectro 20D Plus RS-232C has spectral bandwidth of 5.0 nm. **This instrument has a single cell holder for test tubes and two adapters for 10 x 10 mm cells.** Spectro 20D Plus RS-232C is **lightweight and compact** with a continuous wavelength from 330 to 1000nm. Equipped with a microprocessor, this instrument is ideal for small laboratories and educational institutions. Spectro 20D Plus RS-232C offers ease of use and reliability which can be used in various applications. Spectro 20D Plus RS-232C can be used extensively for qualitative and quantitative analysis.

Spectro 20D Plus RS-232C is a **versatile digital spectrophotometer** with a built-in RS-232C port permits the instrument to be linked to a computer and a printer to display the photometric and spectral data on the PC monitor. This interface provides easy communication between the spectrophotometer and the Window based PC.

Spectro 20D Plus RS-232C utilizes **soft keys for user-friendliness**; with buttons to access settings such as absorption, transmission, concentration and factor, adjusting 0A and 100%T. Spectro 20D Plus RS-232C is able to send data to the computer or print directly from the Spectro to an external printer.

Labomed, Inc. is certified by ISO-9001-2000, has CE Conformity and is FDA Licensed.

Features

- Microprocessor controlled.
- Large Digital LCD Display.
- Auto transmission, absorption, concentration and factor.
- Silicon Photodiode detector
- It can use square or round cells.
- Interface RS232 to link to a computer.
- Wide continuous wavelength range.
- Built in stray light filter.
- 1200 /mm. grating mirror.
- Computer System is optional (NOT INCLUDED).
- Can use 13mm. test tube.

Accessories

- 2 adaptors for 10x10 mm. square cells
- 2 optical glass cells 10x10 mm.
- Dust cover.
- Instruction Manual.
- 1 Zero absorption cell.
- 1 Power Cable

Technical Specifications

Display:	LCD	Absorption Range:	0-2A.
Wavelength Range:	330-1000nm.	Concentration Range:	0-2000
Light Source:	Tungsten Halogen Lamp	Transmittance Range:	0-150% (T)
Detector:	Silicon Photodiode	Photometric Accuracy:	1% T
Wavelength Accuracy :	1nm.	Monochromator:	1200 Lines/Grating Mirror
Wavelength Reliability:	1nm.	Dimensions:	12"(W) x 10" (D) x 8" (H) (inches)
Spectral Bandwidth:	5nm.	Net Weight :	7 kg.